


Waters Alliance Upgrade Kit

2690/2695/2790/2795

The Biotech Systec AF[®] Degassing Replacement Kit for Waters Alliance HPLC instrument is simple to install in half an hour by any laboratory technician.

Please note that Waters degassing Chambers have a slightly different outer configuration but the inside of the chamber is, including the degassing membrane, identical. Biotech has designed a unique bracket which can be mounted into the Alliance Degasser. The unique design of the bracket allows mounting of either Waters original chambers, our Systec AF[®] chambers or a combination of both if there is no need to change all chambers.

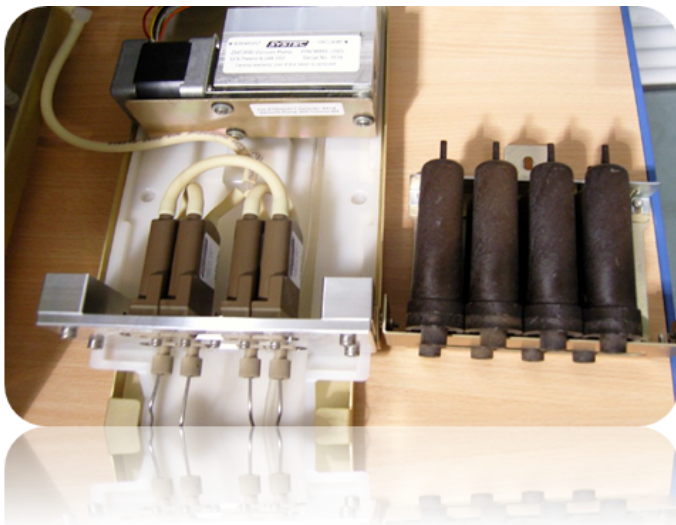


On-line vacuum degassing has been shown to be more effective at removing dissolved gas from mobile phases than other methods. By upgrading, you also take advantage of the latest vacuum degassing technologies.

Contact any of our offices bellow for prices and technical assistance

Biotech Sweden AB
 info@biotech.se
 www.biotech.se
 +46 300 56 91 80

Biotech USA
 bmathieu@biotechusa.us
 www.biotechusa.us
 612-703-5718

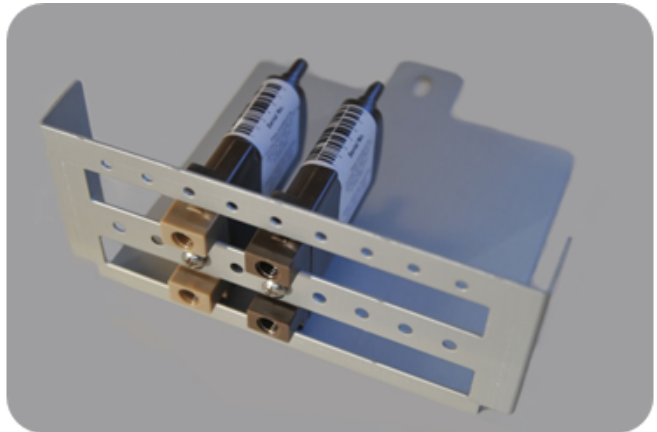

The Waters Degasser Upgrade Kit is complete and consists of the following parts:

- ◆ P/N 18385002 Bracket for Mini Chamber for Waters Alliance HPLC
- ◆ P/N 9000-1004, 480 μ L Systec AF[®] Mini Degassing Chamber
- ◆ P/N 9000-1184, 480 μ L Systec AF[®] Mini Degassing chamber. Recommended for GPC and normal phase applications. (THF and HFIP resistant)
- ◆ P/N 6130-1053. (0,5 M Santoprene Tubing)
- ◆ P/N Bio.1.2. Barbed Connectors for the bracket
- ◆ P/N Bio.1.1. Red plug

